



PROSE Device Use and Care Guide

www.bostonsight.org



BostonSight
464 Hillside Ave, Ste 205
Needham, MA 02494

781-726-7337

TABLE OF CONTENTS

Contact Information..... 3

SET UP and APPLICATION

Set Up: PROSE Device Application..... 4

Preservative-Free Saline Protocol..... 4-5

Key Points for Optimal PROSE Device Application..... 6

The Lid Spread is Everything..... 7

How to Apply Your PROSE Device with a Plunger..... 8-9

How to Apply Your PROSE Device with a Lighted Plunger..... 10

How to Apply Your PROSE Device with a Plunger Stand..... 11

Bubbles Under Your PROSE Device..... 12

REMOVAL

Set Up: PROSE Device Removal..... 13

How to Remove Your PROSE Device with a Plunger..... 13

How to Remove Your PROSE Device without a Plunger..... 14

DISINFECTION and CARE

Set Up: Your PROSE Case..... 15

Set Up: PROSE Device Care and Cleaning..... 16

Nightly Disinfection Protocol for Your PROSE Device..... 17-18

Morning Protocol for the Application of Your PROSE Device..... 19

FAQ for Current BostonSight PROSE Patients..... 20-26

PROSE Solution and Accessory Buying Guide..... 27-31

CONTACT INFORMATION

www.bostonsight.org
464 Hillside Avenue, Suite 205
Needham, MA 02494

Main Phone: (781) 726-7337
Main Fax: (781) 726-7310
Doctors and staff CANNOT respond to medical issues via email.

PROSE device problems including application, removal, or cleaning

1. Review this guide
2. Visit www.bostonsight.org/patient-resources to view our Device Application and Removal video series.
3. Call 781-726-7494 to leave a message for an ophthalmic technician or email triage@bostonsight.org

Patient Emergencies (BostonSight Needham only)
(after hours, weekends, or holidays)

Please note the patient emergency line is for BostonSight Needham patients ONLY. If you are not a BostonSight Needham patient, please contact your PROSE provider directly.

Patient Emergency Line: (617) 755-5929

The emergency telephone line should NOT be used in cases of broken or lost devices, medication questions, or solution and cleaning questions. Please call our office during normal business hours to discuss any non-emergency issues you may be having.

Order a replacement PROSE device(s) or for damaged or lost devices, please contact:

Device Coordinator: (781) 726-7476

If you are requesting device servicing, complete the form linked below and include it with your devices.
<https://www.bostonsight.org/patient-resources>

Insurance Pre-Authorization

Loren Stead: (781) 726-7435

Insurance claim follow-up/appeal

Lynn Carter: (781) 726-7505

Warning

Increased eye redness, pain, discharge, sensitivity to light, and/or significant reduction in vision may signal infection and requires IMMEDIATE attention. Please contact your primary eye care provider and ask them to keep your PROSE doctor aware of any findings.

application

Set Up: PROSE Device Application

You will need a soft towel to be placed on top of a clean work surface where you will apply your device(s).

You will need the following supplies:

- An application plunger
- A removal plunger
- A small mirror
- Preservative-free saline

Preservative-Free Saline Protocol

Filling or Rinsing the Device

Below are 5 options for filling and rinsing PROSE devices:

- ① PuriLens Plus Preservative-Free Saline is a sterile saline solution buffered with boric acid and sodium borate. It is not a stand-alone cleaner or disinfection solution and it does not substitute for nightly cleaning and disinfection. **WARNING:** PuriLens Plus Solution is not to be used with the PuriLens cleaning chamber for PROSE. This chamber is to be used for other contact lenses, but not PROSE devices. PuriLens Plus solution can be ordered directly from the manufacturer at www.purilens.com or by calling (800) 622-0777. It is also available at www.amazon.com.
- ② ScleralFil by Bausch + Lomb is a sterile, buffered isotonic saline solution. Each carton contains 30 10ml vials. It is available at www.amazon.com or www.dryeyeshop.com.
- ③ Nutrifill is a preservative-free buffered saline with natural electrolytes. It comes in a box of 35, 10ml ampules and is available at www.nutrifill.com, www.myeyesupply.com, or www.dryeyeshop.com.
- ④ Sodium Chloride 0.9% vials for nebulization. This is unbuffered, preservative-free saline that comes in 3ml or 5ml vials. An example is the AddiPak Sodium Chloride 0.9% vials. You can request a prescription, however they may not be covered because they are intended for nebulizers and are used off label. They are also available for purchase without a prescription at www.medonthego.com, www.amazon.com, www.myeyesupply.com, or www.dryeyeshop.com.
- ⑤ LacriPure is an unbuffered preservative-free saline, available as 5ml vials in a 98-count box. It is available at www.meniconwebstore.com or www.dryeyeshop.com.



Preservative-Free Saline Protocol, continued

Rinsing Only


To save money, an option for **RINSING ONLY** is to use Bausch + Lomb Sensitive Eyes Plus Saline Solution. This solution contains preservatives and under NO circumstances should be used to fill PROSE devices. When using this solution for rinsing, your PROSE devices should be thoroughly rinsed with preservative-free saline.




application — key points

Key Points for Optimal BostonSight PROSE Device Application

- ① Wash your hands with non-moisturizing soap and dry them completely with a lint-free towel before you begin.
- ② Gently handle only the outside edges of the device.
 - Your PROSE devices may have dots. For patients with devices for both eyes, the RIGHT device is marked with ONE black dot. The LEFT device will have either NO dots or TWO black dots.
 - When either device has one or two black dots, the device should be applied with the dot(s) oriented at the 12 o'clock position, unless otherwise instructed by your doctor.


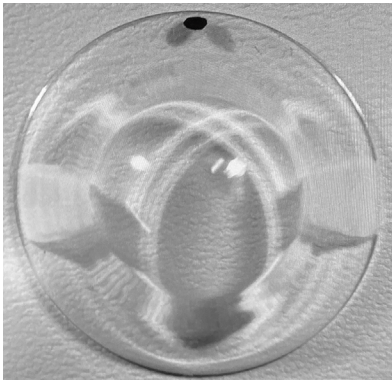


1 dot = Right

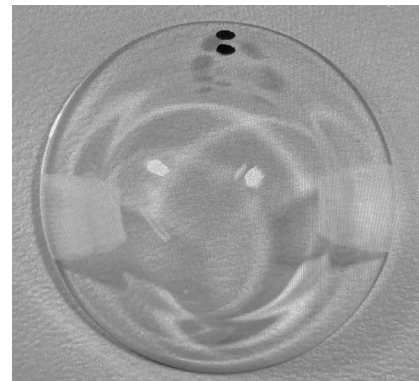

2 dots = Left

- ③ Your eye must be open WIDE – the lid spread is everything! See next page for examples.
- ④ The PROSE device must be centered within the lid spread.
- ⑤ The plunger must be vertical. Your PROSE device must retain the preservative-free saline in the device to avoid bubbles.
- ⑥ Overfill your device to improve your chances of a bubble-free application. Let the saline drip, not stream, into the device bowl.

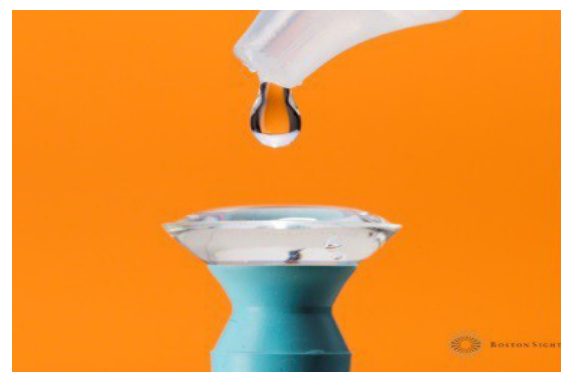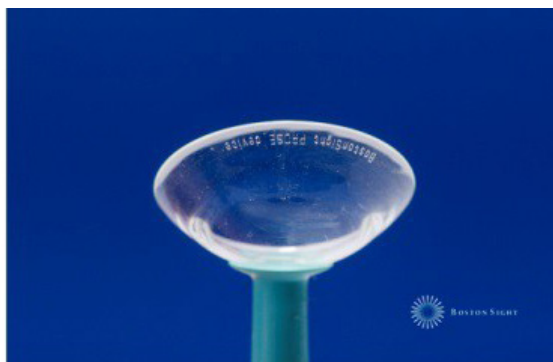

The Lid Spread is Everything



Incorrect! Fingertips are not centered.
Lid spread is not a perfect circle.



Incorrect! Lids slipped from fingertips.
Lids are not open wide enough.




Perfect! A round opening with plenty of
clearance for the device application.

application – plunger

How to Apply Your PROSE Device With a Plunger


- ① Empty the contents of the case (the PROSE device) into the palm of your hand. With the other hand, GENTLY grasp the device by the edges. Note: the device should be wet during this process.
- ② Squeeze in on the sides of the plunger before coming into contact with the device.
- ③ Place the device over the face of the plunger and release side pressure on plunger, creating suction between the plunger and the device.
- ④ Release the sides of the device, holding it with the plunger.

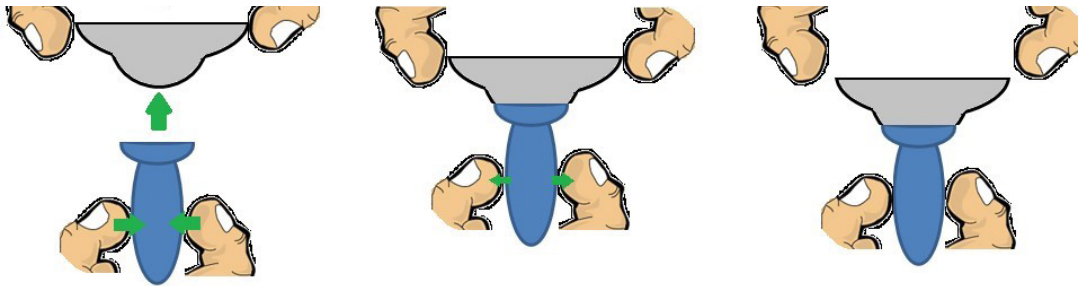

IMPORTANT: Holding the device with the fingers while applying the plunger lessens the force against the device. **ALWAYS** handle the device by the edges when attempting to attach the plunger to the device. **DO NOT** attach the plunger while the device is resting on a hard surface such as a table or counter top. This will place excessive force on the device, risking breakage. Never squeeze the sides or place pressure on the face of the device.

IMPROPER USE OF PLUNGER DURING DEVICE APPLICATION

- DO NOT apply force with plunger in a downward motion causing device to break
- DO NOT use the plunger to contact the device while it is sitting on a hard surface
- DO NOT attempt to take a device out of the case with a plunger

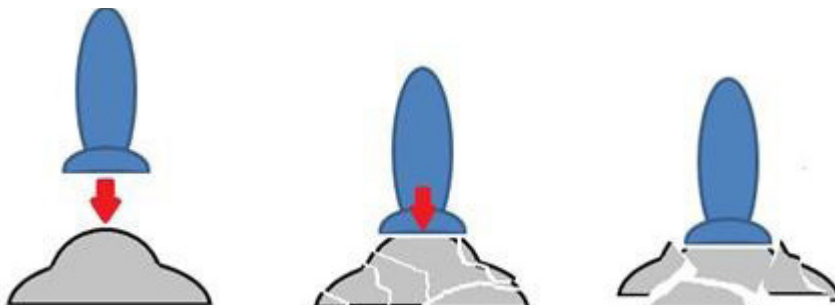

How to Apply Your PROSE Device With a Plunger, continued


- ⑤ Holding the plunger vertically, overfill the device with preservative-free saline. Let the solution drip, not stream into the device, in order to prevent bubbles from forming.
- ⑥ Maintain proper lid spread. Hold the eyelids at the edges, where the lashes meet the lids. The contact points should be centered on both lids to form a perfectly round opening with plenty of clearance for the device.
- ⑦ Tuck your chin to your chest. Look straight down. Keep the plunger centered in your eye by focusing on the black hole in the center of the plunger.
- ⑧ Continue approaching your eye smoothly and steadily, and stare at the black hole as if you are looking through it. Some people find it useful to shut the eye that they are not applying the device to.
- ⑨ You will feel the fluid, but you are not there yet. Press your device onto the eye with firm pressure.
- ⑩ Capture the device with the eyelids by releasing the lids. Don't squint. Maintain the upward pressure on the device while squeezing the plunger in order to release it from the device.
- ⑪ The device has now been released onto your eye! Continue squeezing the plunger while you pull the plunger away from the eye.
- ⑫ Now check your device in the mirror. Shine a penlight or mobile phone light, at an angle, onto the surface of the device and look for the presence of any bubbles. If you detect bubbles the device **MUST** be removed and reapplied.
- ⑬ After application of both devices, clean the tip of your plunger thoroughly with an alcohol prep pad and allow to dry completely. Replace your plungers every 3 months.



application – lighted plunger

How to Apply Your PROSE Device With a Lighted Plunger

- ① Wet the fenestrated end of the application plunger with preservative-free saline. Attach the plunger to the plunger light by inserting the light into the hole at the end of the fenestrated plunger.
- ② Hold a cleaned and disinfected PROSE device gently by the edges. Apply to an application plunger by gently squeezing the plunger at the line while placing the device on the cup and releasing the 'squeeze' to create suction. (Refer to page 8 for complete instructions.) Turn on the green light by pressing the black button.
- ③ Holding the plunger vertically, overfill the device with preservative-free saline. Let the solution drip, not stream into the device, to prevent bubbles from forming.
- ④ Maintain proper lid spread. Hold the eyelids at the edges, where the lashes meet the lids. The contact points should be centered on both lids to form a perfectly round opening with plenty of clearance for the device.
- ⑤ Tuck your chin to your chest. Look straight down. Keep the plunger centered in your eye by focusing on the green light in the center of the plunger.
- ⑥ Continue approaching your eye smoothly and steadily, and stare at the green light as if you are looking through it. Some people find it useful to shut the eye they are not applying the device to.
- ⑦ You will feel the fluid but you are not there yet. Press your device onto the eye with firm pressure.
- ⑧ Capture the device with the eyelids by releasing the lids. Don't squint. Maintain the upward pressure on the device while squeezing the plunger in order to release the plunger from the device.
- ⑨ The device has now been released onto your eye! Continue squeezing the plunger while you pull it away from your eye.
- ⑩ Check your device in the mirror. Shine a penlight or mobile phone light, at an angle, onto the surface of the device and look for the presence of any bubbles. If you detect bubbles the device **MUST** be removed and reapplied.
- ⑪ After application of both devices, clean the tip of your plunger thoroughly with an alcohol prep pad and allow to dry completely. Replace your plungers every 3 months.

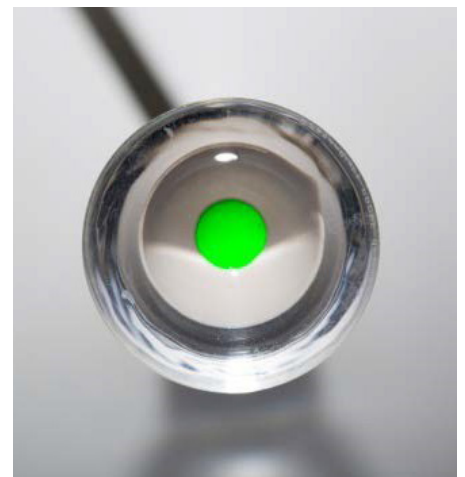

application – plunger stand

How to Apply Your PROSE Device With a Plunger Stand With or Without a Lighted Plunger

- ① If using a light, attach the fenestrated plunger to the plunger light by inserting the light into the hole at the end of the plunger. Next, place the plunger in the hole at the top of the stand and feed the wire through the slot in the stand. Turn on the light. Adjust the height of the stand, if necessary, by placing on top of something such as a book or a tissue box.
- ② Hold the cleaned and disinfected PROSE device gently by the edges, and place onto the application plunger. Do NOT squeeze the plunger to attach the device; this will allow it to release easily when applying.
- ③ Overfill the device with preservative-free saline. Let the solution drip, not stream into the device, in order to prevent bubbles from forming.
- ④ Hold the eyelids at the edges, where the lashes meet the lids. The contact points should be centered on both lids to form a perfectly round opening with plenty of clearance for the device.
- ⑤ Tuck your chin to your chest. Look straight down. Keep the plunger centered in your eye by focusing on the black hole or green light in the center of the plunger.
- ⑥ Begin to approach the device by moving your head downward smoothly and steadily. The approaching device and plunger must be centered within the opening and your focus must remain on the black hole or green light. Some people find it useful to shut the eye that they are not applying the device to.
- ⑦ You will feel the fluid but you are not there yet. Continue the downward motion until the device is pressed firmly against your eye.
- ⑧ Capture the device with the eyelids by releasing the lids.
- ⑨ The device has now been released onto your eye! Keep your eye closed while you sit up, backing away from the plunger with the device on your eye.
- ⑩ Now check your device in the mirror. Shine a penlight or mobile phone light, at an angle, onto the surface of the device and look for the presence of bubbles. If you detect bubbles the device MUST be removed and reapplied.
- ⑪ After application of both devices, clean the tip of your plunger thoroughly with an alcohol prep pad and allow to dry completely. Replace your plungers every 3 months.



application – bubbles

Examples of Bubbles Under Your PROSE Device

Single bubble:



Multiple bubbles:



If you see a bubble, the device **MUST be removed and reapplied.**

Set Up: PROSE Device Removal


You will need to have a soft towel to be placed on top of a clean work surface where you will remove your device(s).

You will need the following supplies:

- A removal plunger
- A small mirror, positioned upright
- Preservative-free saline

How to Remove Your PROSE Device with a Plunger

- ① Position the mirror upright on your work surface.
 - ② Lubricate your eye thoroughly with several drops of preservative-free saline.
 - ③ Wet a removal plunger with the preservative-free saline.
 - ④ Look straight ahead into the mirror.
 - ⑤ Just as when applying your PROSE device, execute proper lid spread. It is important to make sure the upper lid is secure during the removal to prevent scratching your eye.
 - ⑥ While maintaining your lid spread, create suction with the plunger at a 45 degree angle. Aim for the 6 o'clock position on the device, at the point where the colored portion of the eye meets the white.
 - ⑦ Press the removal plunger firmly but gently against the device. Break the suction by tipping the bottom of the plunger up and then pull the plunger away in an arc-down-and-out manner, taking care not to scratch your cornea. On some occasions, you may need to pull and relax the plunger rapidly and repeatedly.
- NOTE: If the plunger is pressed too close to the center of the device, directly over the pupil, it may not be possible to break suction and remove the device.
- ⑧ Maintain the lid spread until the device is out of your eye.
 - ⑨ Hold the device gently by the edges and remove the plunger from the device by twisting and maneuvering it towards the edge of the device.
 - ⑩ Clean the tip of your plunger thoroughly with an alcohol prep pad and allow to dry completely. Replace your plungers every 3 months.


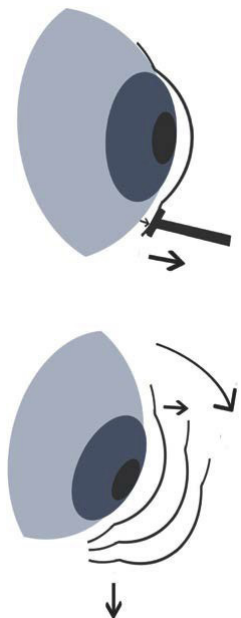


removal – without a plunger


How to Remove Your PROSE Device Without a Plunger

① Lubricate your eye thoroughly with several drops of preservative-free saline.


② Press the tip of an index finger on the very edge of your upper lid at the 12 o'clock position and gently push it into your eye just above the top edge of the device. (If you have deep-set eyes, this will require you to reach under your brow and look slightly down.)

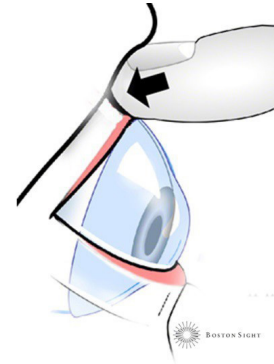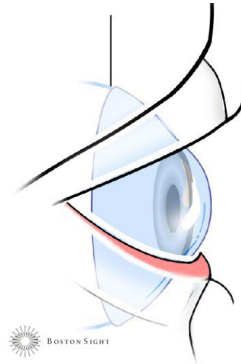


③ Block the device from moving down by pressing the tip of the index finger of your other hand on the very edge of your lower lid just below the bottom edge of the device.

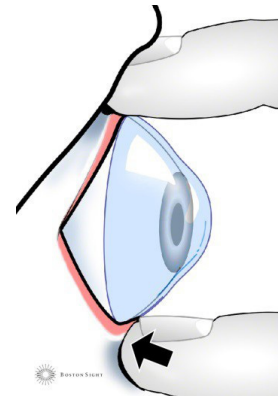

④ Look upward until you feel the device being blocked by pressure on the upper lid margin.

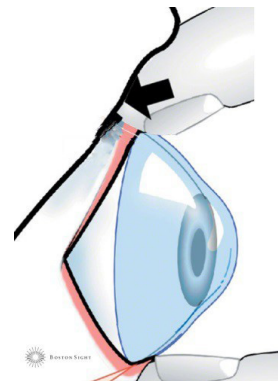

⑤ Scoop the upper edge of the device off the eye with the upper lid margin, look up further, and blink slowly to remove the device.

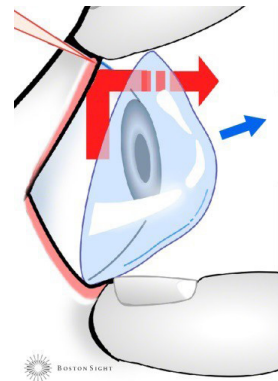

Set Up: Your PROSE Case

To start, you will need to have a clean work surface and the following supplies:


- One PROSE case
- One unopened package of Clear Care
- Please note that you can use original Clear Care or Clear Care Plus



- ① Remove the lid from the PROSE case.
- ② Remove the lid from the small Clear Care case.
- ③ Remove/break off the platinum disc from the arbor on the bottom of the Clear Care case (as illustrated below).
- ④ Discard the Clear Care case.
***NEVER attempt to disinfect your PROSE device in the case included with the Clear Care Solution**
- ⑤ If your PROSE case still has the arbor (stem) attached to the bottom of the basket assembly, it will need to be removed by using pliers to grasp the arbor and gently bending the arbor until it breaks. If not, continue to step 6.
- ⑥ Drop the platinum disc from the Clear Care case into the PROSE case so that it sits on the bottom of the PROSE case.

We recommend replacing the platinum disc every 2 months.

- ⑦ Your PROSE case is now properly set up for nightly disinfection.

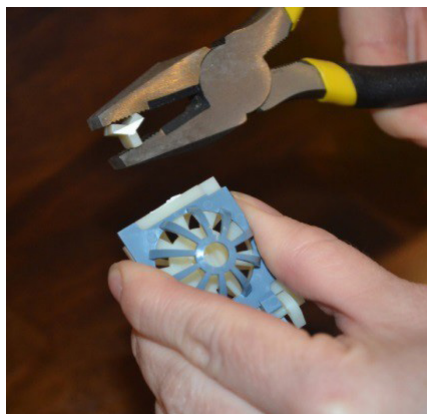

We recommend replacing your PROSE case every 3 months.

disinfection and care

Set Up: PROSE Device Care and Cleaning

You will need a soft towel to be placed on top of a clean work surface where you will clean your device(s).

You will need the following supplies:

- A PROSE case with platinum disc from the Clear Care disinfection case
- Preservative-free saline
- MiraFlow, Sereine, Lens Fresh, Unique pH, Boston Simplus or Tangible Clean* (*Tangible Clean is for Hydra-PEG devices only)
- One bottle of Clear Care or 3% hydrogen peroxide

Nightly Disinfection Protocol for Your PROSE Devices


- ① Wash your hands thoroughly with non-moisturizing soap and dry completely with a lint-free towel before you begin.
- ② Place your device in the palm of your hand. Place one drop of your daily cleaner on the outside surface of the device and rub for 15 seconds. Place another drop on the inner surface and then rub again for 15 seconds.

DO NOT rub the device between your fingers to clean it. This can cause the device to warp.

- ③ While handling the PROSE device gently by the edges, rinse the device thoroughly with either preservative-free saline or Bausch + Lomb Sensitive Eyes Plus. **Under no circumstances should tap water be used for rinsing.**

- ④ Gently separate the basket from the basket assembly by pulling the tab located near the arbor so that the basket is horizontal.

- ⑤ Place the device, bowl side up, into the basket (see photo at right for correct placement). While keeping the basket horizontal, bring the basket assembly down onto the basket. **To avoid damaging or losing your PROSE devices, DO NOT attempt to bring the basket up to the basket assembly.**

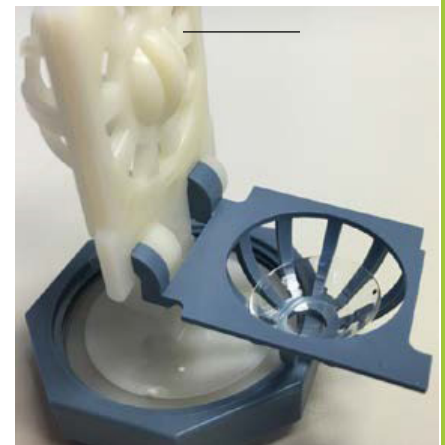


- ⑥ Repeat for the left device and place into the empty basket.

- ⑦ Fill the PROSE case with Clear Care solution or 3% hydrogen peroxide (H_2O_2) to the fill line. It is okay if the platinum disc floats up to the top – the basket will push it back down.

- ⑧ Place the basket assembly in the case and secure the cover. **DO NOT** over tighten the cover because the threads may split and break. It is okay if the platinum disc does not sit flat at the bottom of the case. It will still work properly even if it is wedged on the side of the case.

To properly allow the Clear Care or 3% hydrogen peroxide solution to neutralize, devices must remain in the case for **AT LEAST 6 hours**.

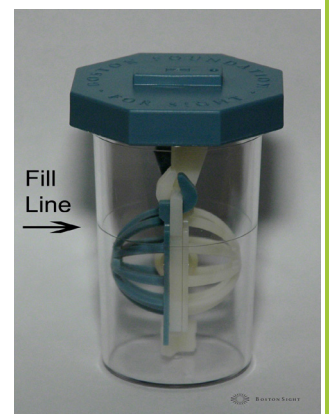

NOTE: If you must remove the devices before the disinfection process has ended, rinse them thoroughly with preservative-free saline. Then soak devices in fresh preservative-free saline, in a clean contact lens case, for at least 5 minutes. Repeat with fresh preservative-free saline solution for an additional 5-minute soak. Rinse again with preservative-free saline before applying to your eye. If using a screw top contact lens case, clean the case with an alcohol wipe and air dry.

disinfection and care p.m.

Nightly Disinfection Protocol for Your PROSE Devices Using 3% hydrogen peroxide (H₂O₂) Without Clear Care Platinum Disc

- ① Clean PROSE device(s) with your daily cleaner by following steps 1-6 on prior page.
- ② Rinse device(s) thoroughly with saline solution.
- ③ Fill the PROSE case with 3% hydrogen peroxide to the fill line.
- ④ Place the basket assembly in the case and secure the cover. **DO NOT** over tighten the cover because the threads may split and break. To properly allow the solution to disinfect, devices must remain in the case for **AT LEAST 6 hours**.
- ⑤ In the morning, remove device(s) from hydrogen peroxide.
- ⑥ Rinse the device(s) thoroughly with preservative-free saline. Then soak device(s) in fresh preservative-free saline, in a clean contact lens case, for at least 5 minutes. Repeat with fresh preservative-free saline solution for an additional 5-minute soak. Rinse again with preservative-free saline before reapplying to your eye.

NOTE: Brown bottle 3% hydrogen peroxide is **NOT FDA** approved for PROSE device disinfection. It is not micro-filtered for impurities and does not contain stabilizing buffers or salts like Clear Care. Some patients report redness and/or discomfort with its use. If at any time while using 3% hydrogen peroxide you experience increased redness, burning, or pain, discontinue usage and use **ONLY** Clear Care.

Morning Protocol for Application of Your PROSE Devices

You will need to have a soft towel to be placed on top of a clean work surface where you will apply your device(s).

You will need the following supplies:

- An application plunger
- A removal plunger
- A small mirror
- Preservative-free saline

- ① After thoroughly washing and drying your hands with non-moisturizing soap and water, remove your device from the PROSE case.
- ② While holding a cleaned and disinfected device gently by the edges, briefly rinse with preservative-free saline. **Under no circumstances should tap water be used for rinsing.**
- ③ Apply device(s) to the eye(s) using your preferred application technique.
- ④ After each use, empty the PROSE case (being careful not to lose the platinum disc) and leave it open to air dry. NEVER clean the PROSE case with alcohol.



NOTE: If you experience a slight stinging sensation upon application of the device in the morning, remove your device from your eye immediately and rinse your eye thoroughly with preservative-free saline. Do this until the stinging sensation dissipates. Rinse your device(s) thoroughly with preservative-free saline and then soak them in fresh preservative-free saline solution for at least 5 minutes before reapplying to your eye. When this occurs, it is likely an indication that the Clear Care platinum disc needs to be replaced.

WARNING: Increased eye redness, pain, discharge, sensitivity to light, and/or significant reduction in vision may signal infection and requires IMMEDIATE attention. Please contact your primary eye care provider and ask them to keep your PROSE doctor aware of any findings.

faq – PROSE device use

Q: What are the daily essentials of my PROSE device(s)? Is there a smaller carrying case available than the one provided?

Daily Essentials include:

- Preservative-free saline
- One clean contact lens case
- One application plunger
- One removal plunger

A small digital camera case or makeup bag, available at any major department store, should fit your daily essentials.

Q: What kind of eye makeup can I use while wearing my BostonSight PROSE device(s)?

Please visit our website for a list of makeup options that are complimentary to PROSE device wear at

www.bostonsight.org/wp-content/uploads/2020/09/Compatible-Makeup-List.pdf

Q: I have two PROSE devices. How can I distinguish between the left and the right device?

For patients with PROSE devices for both eyes, the RIGHT device is marked with ONE black dot. The left device will have either NO dots or TWO black dots. The device should be applied with the dot(s) oriented at the 12 o'clock position, unless instructed otherwise by your doctor.

Q: Can I sleep or nap in my PROSE device(s)?

You can nap for a maximum of 30 minutes. **DO NOT** sleep with your PROSE devices, as this can greatly increase the risk of an eye infection.

Q: What should I do with my PROSE device(s) if I remove them for a short period of time to shower, swim, or nap?

Devices should be placed in a screw top contact lens case with either Unique pH or Boston Simplus. Before reapplying, the PROSE device(s) should be cleaned with your daily cleaner and rinsed well with preservative-free saline. The case should then be emptied, wiped out with an alcohol prep pad, and allowed to air dry.

Q: Can I swim or shower while wearing the PROSE device(s)?

No, you cannot swim or shower while wearing PROSE devices. Exposure to impure and non-sterile water (tap, lake, pool, etc) while wearing your PROSE device increases your risk of a sight-threatening infection. For any further specific questions regarding this, please speak to your PROSE doctor.

faq – PROSE device use

Q: I just received my PROSE device(s) in the mail. How do I prepare them for application?
PROSE devices will be shipped wet in a case. Upon receiving your device(s), please follow the steps below:

1. Wash your hands thoroughly with non-moisturizing soap and water before handling your PROSE device(s).
2. Clean your PROSE device(s) as instructed below:
 - Place one PROSE device in the palm of your hand.
 - Place one or two drops of your daily cleaner on the outside surface of the PROSE device and rub for 15 seconds. Flip the PROSE device over, place another one or two drops of your daily cleaner, and rub the inside surface for 15 seconds.
3. While handling your PROSE device gently by the edges, rinse with the sterile, preservative-free saline.
4. Your PROSE device is ready for application.

Q: My condition causes mucus to accumulate on the outside of my PROSE device(s). What can I do?

Under no circumstances should tap water be used for rinsing.



For any **plasma-treated** PROSE device, you can use a few drops of sterile preservative-free saline on the end of a Q-tip to “squeegee” the outside surface of the device. Alternatively, PROSE devices can be refreshed: remove, clean with your daily cleaner, rinse thoroughly with preservative-free saline, and reapply.

Hydra-PEG treated devices: Any Hydra-PEG treated PROSE device can be refreshed: Remove, clean with Hydra-PEG approved daily cleaner, rinse thoroughly with preservative-free saline, and reapply. Do not “squeegee” a Hydra-PEG device since this could damage the integrity of the coating.

Q: I use medicated eye drops and wear PROSE devices. How do I use the drops?

If the medication is used once or twice per day, put the eye drop in before PROSE device application and after device removal, making sure to wait 10-15 minutes prior to application of PROSE device. If the medication is three times per day or more, then consult with your primary eye care provider to determine whether it is acceptable to apply drops while wearing the PROSE device, as the absorption is compromised.

faq – traveling with PROSE devices

Q: I am traveling to _____ country/non-US city. Can you please tell me where there is a pharmacy where I can buy the solutions I need for my PROSE device(s)?
Unfortunately, we have no way of knowing which solutions are available in which countries. We recommend that for extended trips, you either pack the solutions that you will need and ship them to yourself at your destination or order some supplies from a U.S. distributor for delivery to your destination address.

Q: What solutions can I take on an airplane?

You can carry on all necessary solutions except 3% hydrogen peroxide and Clear Care, which are not allowed on an aircraft under any circumstances. These must be packed in your checked luggage. We will provide you with a medical necessity letter addressed to the Transportation Security Administration (TSA) detailing what you should be allowed to bring onto the airplane due to medical reasons.

Q: How often do I need to be seen for follow-up care by my BostonSight PROSE provider?

Once the 6-month global treatment period is complete, patients are expected to return annually to their BostonSight PROSE provider for a comprehensive PROSE treatment evaluation, or sooner if instructed by your doctor.

Q: Do I continue to see my regular ophthalmologist or optometrist?

Yes. You should follow up with your referring eye doctor within 3 to 4 weeks after completing PROSE treatment. Continue to see your primary eye care provider for complete eye exams, including dilated examinations, on an annual basis or per their recommendation. BostonSight does not perform routine eye care and vision exams.

Q: What should I do if I experience any of the following symptoms?

Discontinue use of your PROSE device(s) and seek medical consultation by your primary eye care provider or emergency facility if you experience any of the following:

- Persistent redness/soreness, never observed before
- Increased sensitivity to light while wearing the device(s)
- Eye pain
- Presence of mucous discharge
- Burning

faq – PROSE device cleaning

Q: What is the most current device cleaning protocol?

We are dedicated to selecting the best products and cleaning regimen for our patients. This means that from time to time we change the products we use. The current cleaning regimen includes:

For Daily Cleaning

- For **plasma-treated** devices: MiraFlow Extra-Strength Daily Cleaner, Sereine Extra Strength Daily Cleaner, or Lens Fresh Extra Strength Contact Lens Cleaner
- For **Hydra-PEG** treated devices: Unique pH, Boston Simplus, or Tangible Clean

For Rinsing

- Preservative-free saline

For Disinfection

- Clear Care or 3% hydrogen peroxide (H₂O₂)

Q: Is H₂O₂ (brown bottle 3% hydrogen peroxide) approved for use with PROSE devices and if so, how do I use it?

Yes, you can use H₂O₂, but its usage comes with a disclaimer.

PLEASE READ ALL OF THE INSTRUCTIONS BELOW BEFORE USING H₂O₂.

Brown bottle 3% hydrogen peroxide is NOT FDA approved for PROSE device disinfection. It is not micro-filtered for impurities and does not contain stabilizing buffers or salts like Clear Care. Some patients report redness and/or discomfort with its use. If at any time while using brown bottle H₂O₂ you experience increased redness, burning, or pain, discontinue usage and use ONLY Clear Care.

The platinum disc that comes packaged with Clear Care (which is not available to purchase by itself) is required to properly disinfect your PROSE devices in the case. Therefore, you must purchase Clear Care at least every 8 weeks to replace the platinum disc.

***NOTE:** At this time, we cannot recommend a generic equivalent of Clear Care available at some retailers (i.e. store brand may say “compares to Clear Care” on the package). Please refer to the BostonSight PROSE Accessory and Solution Buying Guide on pages 27-30.

faq – PROSE device cleaning

Q: Is there a cheaper alternative to the preservative-free saline to fill my device?

Yes. You may use Sodium Chloride Inhalation Solution USP, 0.9% as a substitute. This is available at your local pharmacy with a prescription or online without a prescription at www.medonthego.com. It is available in 3, 5, and 10 mL vials and comes in boxes of 100–144.

NOTE: While most people can use Sodium Chloride Inhalation Solution USP, 0.9% in place of the preservative-free saline, some will develop mild redness, burning, and/or discomfort, as it is not buffered. If you develop these symptoms, discontinue use immediately and resume use of a buffered preservative-free saline. **YOU MUST always FILL device(s) with preservative-free sterile saline solution.**

Q: Is okay to use the drops that claim to reduce protein buildup on contact lenses on my PROSE device(s)? If so, which brand should I use?

****Enzymatic cleaner is NOT approved for devices treated with Hydra-PEG.**

We suggest using Ultrazyme Enzymatic Cleaner or Boston One Step Liquid Enzymatic Cleaner as directed on the package inserts.

- If using Ultrazyme, add one tablet to your PROSE case with Clear Care or 3% hydrogen peroxide solution and disinfect normally as instructed in this guide.
- If using Boston One-Step, empty one vial into your PROSE case with Clear Care or 3% hydrogen peroxide solution and disinfect normally as instructed in this guide.

Q: How often do I need to replace my PROSE case and plungers?

With proper use, the plungers should be replaced every 3 months; sooner if the edges become rough, uneven, or the suction is no longer sufficient. The PROSE case should be replaced every 3 months. PROSE cases and plungers are available at The Dry Eye Shop: www.dryeyeshop.com/

Q: Where can I purchase my cleaning supplies?

Please use the link to the comprehensive BostonSight PROSE Solution and Accessory Buying Guide for more information. This guide can be found on pages 27-30.

faq – PROSE device care

Q: How do I store my backup PROSE device(s)?

For plasma-treated PROSE device WITHOUT Hydra-PEG coating, the device can be stored per the manufacturer's recommendation in your PROSE case untampered with (no opening/closing) if not worn immediately after overnight disinfection.

When ready to wear, clean with Miraflow, Lens Fresh, Sereine, Unique pH or Boston Simplus. Rinse the device thoroughly with preservative-free saline before applying. If storing for more than 1 week, we recommend removing the device from the PROSE case and drying with a soft cloth. When completely dry, the PROSE device can then be stored DRY (without solution) indefinitely in a clean screw top contact lens case. When you want to use the PROSE device again, make sure to disinfect overnight in Clear Care prior to use.

For Hydra-PEG treated PROSE devices, be sure to store wet if not worn after overnight disinfection. Devices can be stored in Unique pH per the manufacturer's recommendation for up to 30 days, if untampered with (no opening/closing). Disinfect overnight in Clear Care prior to wearing your device(s).


Q: What is involved in getting my PROSE devices serviced? (Servicing includes plasma treatment or reapplication of Hydra-PEG and darkening of dots.)

Your PROSE devices are plasma treated or have Hydra-PEG coating applied as the last step of their manufacturing. The treatments make them more wettable so that vision and comfort are both improved. Simple wearing and handling, as well as the care solutions, can wear down the treated surface. This can make the PROSE device(s) cloudy, requiring more frequent blinking, rewetting and 'squeegeeing' to clear the vision. They can also become increasingly more uncomfortable to wear because of increased dryness. As to how often they may need treatment, this varies greatly, from not at all to every 3 to 4 months.

If you suspect or have been told that your PROSE device(s) need to be plasma treated:

- Complete a PROSE Device Service Form at www.bostonsight.org/patient-resources/
- Include a check for payment or pay at www.bostonsight.org/patient-resources/
- **Send them to us 'Attention: Device Coordinator'**
- We recommend you send PROSE device(s) via FedEx, UPS, or by other means that allows tracking. This will reduce the chances that they may be lost. Once we receive them, we will treat them and send them back to you.
- If you choose to have your PROSE device(s) returned to you via USPS regular mail, there is no associated shipping cost. Requests for expedited shipping using FedEx will incur a charge. Please note that payment for FedEx shipping must be made by credit card before your device(s) are returned. Without exception, all packages, regardless of shipping method, require a signature on arrival.

PROSE accessory buying guide

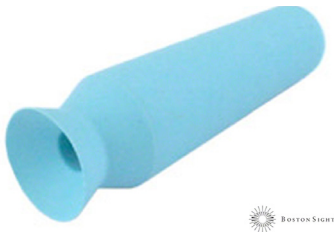

Application Plunger (DMV Scleral Cup)
www.dryeyeshop.com
877-693-7939



See-Green® Plunger (Fenestrated Application Plunger)
www.seegreenlight.com/
413-563-2477



Removal Plunger (DMV Ultra Contact Lens Remover)
www.dryeyeshop.com
877-693-7939



See-Green® Scleral Lens Inserter (Light)
www.seegreenlight.com
413-563-2477
*register your See-Green Lens Inserter for a 2-year manufacturer's warranty at www.seegreenlight.com/registration



See-Green® Scleral System (Light and Stand)
www.seegreenlight.com
413-563-2477



PROSE Case
www.dryeyeshop.com
877-693-7939

Both Dry Eye Shop and See-Green ship internationally.

PROSE solutions

Daily Cleaning Solutions Plasma-treated lenses



Lens Fresh



MiraFlow



Sereine

Daily Cleaning Solutions Hydra-PEG treated lenses



Boston Simplus



Tangible Clean



Unique pH


Rinsing and Filling Solutions (preservative-free salines and alternatives)



PuriLens



LacriPure

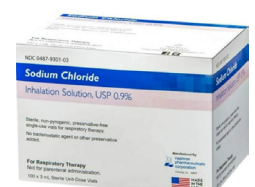

Sodium Chloride 0.9%
for Inhalation



ScleralFil



Nutrifill

Rinsing Only Solutions



B+L Sensitive Eyes

Disinfecting Solutions



Clear Care



3% hydrogen peroxide

Solutions buying guide

Bausch + Lomb Sensitive Eyes Plus Saline

Available at the following retail locations: CVS, Walgreens, Target, Walmart, and most grocery stores and at www.amazon.com, www.cvs.com, www.walgreens.com, www.target.com, www.theonlinedrugstore.com, and www.walmart.com

BostonSimplus

Available at the following retail locations: CVS, Walgreens, Target, Walmart, and most grocery stores and at www.amazon.com, www.cvs.com, www.walgreens.com, www.target.com, www.theonlinedrugstore.com, www.walmart.com and www.dryeyeshop.com

Clear Care

Available at the following retail locations: CVS; Walgreens; Rite Aid; Target, Walmart, and most grocery stores and online at www.cvs.com, www.walgreens.com, www.walmart.com, www.theonlinedrugstore.com, www.amazon.com, and www.aclens.com

LacriPure

Available at www.meniconwebstore.com and www.dryeyeshop.com

Lens Fresh Extra Strength Contact Lens Cleaner

Available at <https://lensfresh.com>

MiraFlow Extra-Strength Daily Cleaner

Available at <http://miraflow.com/purchase-here>

Nutrifill

Available at www.dryeyeshop.com and www.nutrifill.com

PuriLens Plus

Available at some Harmon Discount Health & Beauty retail locations and online at www.purilens.com, www.amazon.com, www.walmart.com, www.aclens.com, and www.dryeyeshop.com

ScleralFil

Available at www.dryeyeshop.com or www.amazon.com

Sereine Extra Strength Cleaner

Available at <https://myeyesupply.com/product/sereine-extra-strength-daily-cleaner/>

Sodium Chloride 0.9% for Inhalation

Available at your local pharmacy with prescription and at www.amazon.com, www.walmart.com, www.medontheho.com, www.myeyesupply.com, and www.dryeyeshop.com

Tangible Clean

Available at www.tangiblescience.com/product/tangible-clean/ and www.dryeyeshop.com

Unique pH

Available at Walgreens and Walmart and at www.amazon.com, www.walgreens.com, www.walmart.com and www.dryeyeshop.com

Buying guide by retailer

www.amazon.com

- Bausch + Lomb Sensitive Eyes Plus Saline
- Clear Care
- Boston Simplus
- Purilens Plus
- ScleralFil
- Sodium Chloride 0.9% for Inhalation
- Unique pH

CVS retail locations and online at www.cvs.com

- Clear Care
- Bausch + Lomb Sensitive Eyes Plus Saline
- Sodium Chloride 0.9% for Inhalation (Rx only)
- Boston Simplus

www.dryeyeshop.com

- Boston Simplus
- Sodium Chloride 0.9% for Inhalation Application Plunger
- LacriPure
- Nutrifill
- PROSE Case
- PuriLens Plus
- Removal Plunger
- ScleralFil
- Tangilble Clean
- Unique pH

www.medonthe-go.com

- Addipak 0.9% sodium chloride for inhalation

www.nutrifill.com

- Nutrifill

www.theonlinedrugstore.com

- Clear Care
- Bausch + Lomb Sensitive Eyes Plus Saline
- Boston Simplus

www.purilens.com

- Purilens Plus

Buying guide by retailer

www.seegreenlight.com

- Fenestrated Application Plunger (See-Green Plunger)
- Lighted Plunger (See-Green Lens Insertter)
- Plunger Stand (See-Green Scleral System)

Target retail locations and online at www.target.com

- Clear Care
- Bausch + Lomb Sensitive Eyes Plus Saline
- Boston Simplus

Walgreens retail locations and online at www.walgreens.com

- Clear Care
- Bausch + Lomb Sensitive Eyes Plus Saline
- Sodium Chloride 0.9% for Inhalation (Rx only, not available online)
- Unique pH
- Boston Simplus

Walmart retail locations and online at www.walmart.com

- Clear Care
- Bausch + Lomb Sensitive Eyes Plus Saline
- Sodium Chloride 0.9% for Inhalation (Rx required for store purchase)
- Unique pH
- Boston Simplus

The following retailers will ship internationally.

- www.amazon.com
- www.seegreenlight.com
- www.dryeyeshop.com
- www.medonthego.com
- www.purilens.com
- www.theonlinedrugstore.com



BostonSight
464 Hillside Avenue, Suite 205
Needham, MA 02494
(781) 726-7337
<http://www.bostonsight.org>